Manual Utilizare

User manual

**RO** Cântar persoane  
MY4806

**GB** Body scale  
MY4806

## Manual de utilizare

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## INSTRUCTION MANUAL

Thank you for buying our bathroom scale, please read this manual carefully before using and maintain it for future reference.

### Product Specification

1. Capacity: 150 kg
2. Graduation: 0.1 kg


### The Battery

The product uses 2 AAA


### To Measure Weight

Before your first use, please follow below steps:

1. Place the scale on a flat horizontal surface. This condition is mandatory – it allows to receive an exact result. Then directly stand on the scale with your feet to turn it on.

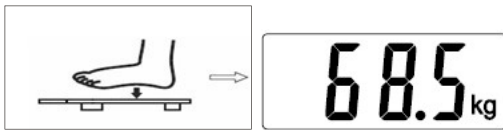2. The factory default unit is “kg”. If not applicable, please set switch on the back of the scale to choose the needed unit “lb” or “st”.

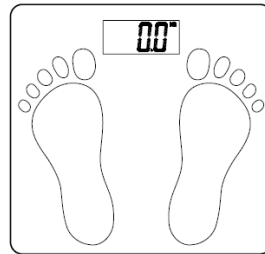3. To stand on the scale directly.
4. The LCD displays your weight, as shown in the below picture.

5. The scale will turn off automatically after 5 seconds without any operations.

### Note

1. “C” will appear if the scale is in unstable conditions. It means it's not the real weight, you need to weight again.

2. “ERR” will appear if the scale is overloaded (Over 150 kg/333LB).

A rectangular digital display showing the text "Err" in a black, monospaced font.

3. “Lo” will appear if the battery power is low and needs to be replaced.

A rectangular digital display showing the text "Lo" in a black, monospaced font.

### Caution

1. Keep the scale off water.
  2. Keep the scale with care.
  3. Do not put it on uneven floor while weighting.
  4. Do not step on the edge of the platform.
  5. Do not step on the scale with wet feet.
  6. Do not use the scale under electromagnetic interference condition.
  7. Please take the batteries out if the scale is not used for a long time.
  8. Do not attempt to open the scale or to remove any components.
- Servicing should be done by qualified technicians only.

### Important

1. The scale you buy is a high-precision electronic instrument. It contains precision sensor, please keep the scale off water and don't drop or jar the scale.
2. Please replace the batteries in time if the device shows low battery. Don't use the new battery mixed with the old ones.

Correct disposal of this Product

(Waste Electrical & Electronic Equipment (WEEE))

Your product is designed and manufactured with high quality materials and components which can be recycled and reused.

This symbol means that electrical and electronic equipment, at their end-of-life should be disposed of separately from your household waste.

Please dispose of this equipment at your local community waste collection/recycling center.

This product complies with European safety and electrical interference directives

This product does not contain environmentally hazardous materials (lead, mercury, cadmium, hexavalent chromium and brominated flammable agents: PBB and PBDE).

CE  