Manual de utilizare

Cântar inteligent

MY4836

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MY4536

Dear users:

Thank you for choosing our Smart Body Analysis Scale.

This product is designed by your body relationship, to monitor the fat, water, muscle, bone, BMI, visceral fat and so on. The product use 2.4G HZ wireless transmission protocol, can connect Bluetooth from mobile or tablet to transmit data.

In order to ensure the accuracy of measurement data and service life, please read the instructions carefully before use and keep it safely for future reference.

Product Features

- 1.Measure Weight, Body Fat, Hydration, Visceral Fat, Muscle Mass, Bone Mass, Calories , BMI and so on.
- 2.Records data to App allowing for you to check progress and maintain goals.

More Details

- Weighting scope/dividing value:5kg-180kg/0.1kg, division compensation:0.05kg
- Battery:4XAAA battery
- Three units: kg/st/lb
- The battery voltage is less than 2.4V,shows “Lo”,15seconds after auto power-off
- Overload display “oL”, Over Value:182.5kg

How to Ensure Accurate Readings

- Always measure barefoot. Do not step on the scale with wet feet, wearing socks or shoes to prevent slipping, damage and inaccurate readings.
- Recommended to be placed on hard and flat surface only.
- Always weigh yourself on the same scale, same place and same surface if possible.
- Recommended measurement time should be after waking up but before the breakfast, and next weight should be two hours after your meal but before your next meal, or two hours after dinner but before bed.
- Erroneous measurements can occur if feet are dirty, dry or extremely callous.
- Body fat percentage estimates will vary with the amount of water in the body, and can be affected by dehydration or over-hydration due to such factors as alcohol consumption, menstruation, illness, intense exercise.
- Stand upright and not to the edge of the scale to ensure good contact with the film or electrodes (the film and electrodes both are conductive media depending scale models).
- Horizontal line should be running below arch of each foot.
- Full feature of scale only available with corresponding app.

Maintenance

- 1.Keep scale surface dry to avoid accident by slipping..
- 2.Please put the scale in a dry place to protect the internal electronic components.

3. Please use micro cloth to clean the scale, don't immersed in water.
4. Please don't remove the scale except replace the battery.
5. Please don't use heavy percussion or falls from a height.

How to use

① Installing the battery

- 1.1 Remove the battery cover at the back of Smart Body Analysis Scale, and insert 4x1.5V AAA batteries (Not included) into the unit. Make sure the polarity is correct when done.
- 1.2 To set "units", press the button located on the back of the scale. Choose the right unit you need, kg , st or st:lb. Then, put back the battery cover.

② Installing App

- 2.1 Switch on Bluetooth
- 2.2 Phone or Tablet connects to Net(WiFi,4G or 3G)
- 2.3 Use your device to scan QR code or search "OKOK International" from App Store or Google Play to install App.

Click "Turn All Categories On" and "Allow" in installing, after some pages click "Start"

③ Pair the devices and register , complete user information

Open "OKOK International"

Enter your email and password, and Click "Join Us", Then you will be asked to set your personal profile, such as nickname, age, gender, and height etc for complete registration.


Completing a personal profile is essential because the data you enter will be used as reference for APP to analyze and assess health situation.

If you have already registered before, you can login to your account by clicking 'Login'.

If you are registered and forgot the Password, please click Login, and select "Forget password". Type your Email, new password and Click Obtain Vcode (abbreviation for verify code). A verification code will be sent to your email. Type the code obtained from your e-mail. Click "Submit" and restore your ID.

④ Pair your body fat scale

After registration you can bind your body fat scale for measurement late. Click "Device" at the bottom right corner, you see the "Binding device screen" illustrated below.

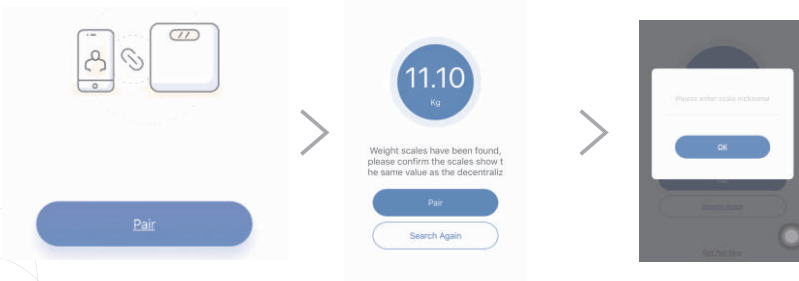
Click "Pair" , step the scale lightly, the number will be displayed on your phone Then click "Pair" ,type a nickname of your body fat scale . After typing nickname, you can step on the scale for measurement.

⑤ Basic Operations

Family Member Management

Enter main screen, click the avatar icon at the top left corner, now all family members' nicknames are display. You can select one member for current measurement results, show or manage his or her data. Also, you can use "Add" operation to add a family member or even a "visitor" for a guest to try a body fat scale. The Member Profile setting screen is shown . You can set personal profile, including nickname, gender, birthday, and height to register to your account.

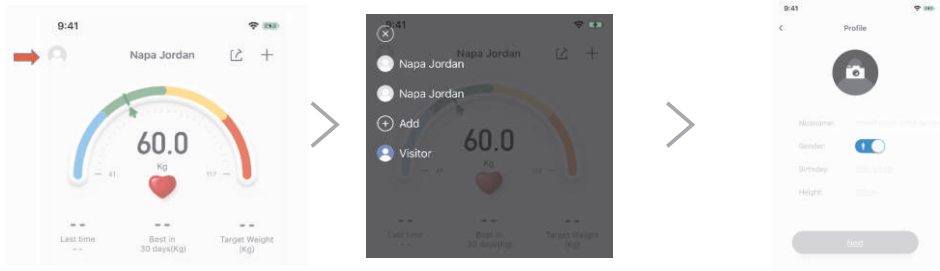
Weight Record Adding Manually

If you like, or you have a usual weight scale rather than a Bluetooth body fat scale, you can add measurement result manually, and use the APP to manage your data. At the main screen, there is a manual adding "+" icon at the top right corner


While you click the "+" icon, a data entering screen is displayed. You can slide the dial to enter the weight, and click "Ok" at the bottom to confirm the number.

Please note, as you only enter weight data, there are NOT other data, such as Fat, Muscle, Visceral Fat, and water etc. The display may look like that shown

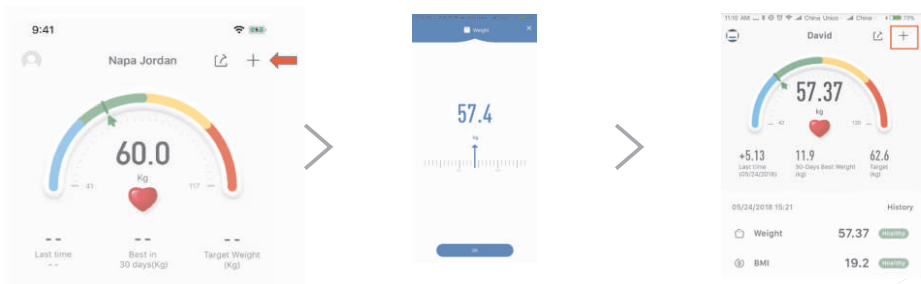Detailed Data Checking

Click on one selected parameter, such as “BMI” on main screen, APP will display detail about BMI

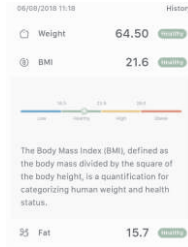Historical Data Checking

Click on “History” on main screen , APP will display historical data of selected family member.

User can slide up or down to check the list of historical records. User can also delete any record. For Android phone, long press one record can delete the record; for IOS phone, slide the record from right to left, a deleting screen is displayed , Click the delete icon to remove the record.

Trend Checking

Click on “Trend” at bottom left corner of main screen, APP will display the diagram of trend of selected parameter, such as weight, BMI, fat, muscle, water, bone mass, or metabolism.


Application Setting

Click on the “Me” icon at the bottom right corner of main screen, then APP will display the “Me” screen with multiply setting operations

Settings Screen

Click the “Settings” button on the Me screen ,App will display the “Settings” screen shown for multiple basic settings.

Here you can do Account Setting, Unit Setting, Weigh Remind and so on.

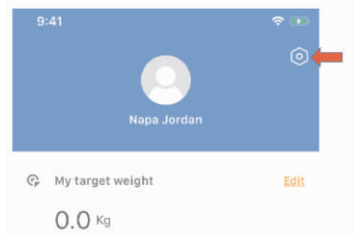More features waiting for you to experience!

Features:

App available for Android & IOS operating systems

Measures weight, body fat, hydration, visceral fat, muscle

mass, bone mass, body mass index (BMI), calories (Kcal)

Conductive surface from indium tin oxide

Temperature and battery level display

Auto on and Auto off

Low battery and overload indicator

Technical Specifications:

Maximum capacity: 180kg

Division: 100 g

Glass size (mm): 300x300x6

LCD size (mm): 59x59

Product size (mm): 300x300x31,5

Battery: 4xAAA (not included)

WARNING: This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

**Correct disposal of this product****(Waste Electrical & Electronic Equipment, WEEE)**

Your product is designed and manufactured with high quality materials and components which can be recycled and reused.

This symbol means that electrical and electronic equipment, at their end of-life should be disposed of separately from your household waste.

Please dispose of this equipment at your local community waste collection/recycling center.

This product complies with European safety and electrical interference directives

This product does not contain environmentally hazardous materials (lead, mercury, cadmium, hexavalent chromium and brominated flammable agents: PBB and PBDE).

Producător:

S.C. COMPLET ELECTRO SERV S.A.

Voluntari- Șos. București Nord, 10

Global City Business Park, Corp 011, et. 10

Telefon relatii clienți: 0374168303

e-mail: relatiiclienti@ces.nt.ro

România